beurer

FS 50

- | | |
|--|---|
| (D) Gesichtssauna
Gebrauchsanleitung | (RUS) Сауна для лица
Инструкция по применению |
| (GB) Facial sauna
Instruction for Use | (PL) Sauna dla twarzy
Instrukcja obsługi |
| (F) Sauna facial
Mode d'emploi | (NL) Gezichts sauna
Gebruiksaanwijzing |
| (E) Sauna facial
Instrucciones para el uso | (P) Sauna facial
Instruções de utilização |
| (I) Sauna per il viso
Istruzioni per l'uso | (GR) Σάουνα προσώπου
Οδηγίες χρήσεως |
| (TR) Yüz saunası
Kullanma Talimatı | |

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Dear customer,

Thank you for choosing one of our products. Our name stands for high quality, thoroughly tested products for the applications in the areas of heat, weight, blood pressure, body temperature, pulse, gentle therapy, massage and air.

Please read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

Best regards,
Your Beurer Team

WARNING

- The unit is only intended for domestic/private use, not for commercial use.
- This device may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely, and are fully aware of the consequent risks of use.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- The mains connection cable of this device cannot be replaced. The device must be scrapped in the event of damage to the cable.
- Check before each use that the device and mains cable are intact. If they are damaged, the device must not be operated.
- Always ensure that the device is securely positioned on a level surface.
Risk of scalding through escaping liquids.
- Warning – the escaping steam is very hot!

1. Signs and symbols

The following symbols are used on the type plate.

-  The device has double protective insulation and therefore complies with protection class 2.

2. Important notes prior to use

- The device is intended only for the purpose described in the instructions for use. The manufacturer cannot be held responsible for damage arising from incorrect or careless use. This device is not suitable for commercial or medical use, but is only intended for private domestic use.
- Repairs must be carried out only by authorised service agents as special tools and special parts are required for repair. Otherwise the guarantee is void.
- Do not use the device if you suffer from heart disease, skin flushing, inflammation, swelling, sensitive areas, diabetes, wounds, burns, eczema or psoriasis.

- During pregnancy or if you have health problems, consult your doctor before using the device.
- Children do not recognise the dangers that can arise in conjunction with electrical devices. Ensure that the device cannot be used by children unsupervised.
- Check before you plug the device into the mains that the local voltage is the same as the voltage given on the device.
- Do not use the device in the vicinity of water (e.g. bath tub, wash basin, etc.).
- Do not use the device in the bath or shower or when you are asleep.
- Do not put it in water or other liquids. Keep the device in a safe place.
- Do not place the device on water-sensitive surfaces.
- Use the device only with the facial sauna steam mask attached.
- Keep the device and mains cable away from hot surfaces.

3. Important instructions during use


- Switch off the device immediately if a fault occurs during use and unplug it from the socket.
- Do not keep your face too close to the mask as otherwise burns can occur because of the hot steam. Close your eyes during use and protect heat-sensitive areas.
- Do not use the device for longer than 20 minutes.
- Stop using it with oils and herbs if skin irritation or anything else unusual occurs in the skin.
- You should also stop using it if you develop pain or feel unwell and you should consult your doctor.
- Do not use the device under covers, towels etc.
- Do not change the accessories while in use as there is a risk of burns if hot liquids spill.
- **Caution:** The overheating protection switches the device off after approx. 2 minute if the water in the water container has run dry.
The device can be started again as normal after a cooling down period of approx. 4 minutes.
- If the device falls into water by accident, unplug it immediately. Have the device checked by a qualified electrician before using it again.

4. Important instructions for after use

- After use, switch off the device and unplug it.
- Never unplug it with wet hands or by pulling on the cable.
- Empty any remaining liquid from the side of the vaporizer. Ensure that no liquids can get into the device through the switch, the cable or the bottom of the device.
- The device should be cleaned after use to avoid the collection of grease and other residues. This ensures optimal operation and prolongs the working life of the device.
- Let the device cool down fully before operating it again.
- If you have any further questions about using our devices, please consult your dealer or the Customer Service department.

5. Device description

1. Facial sauna steam mask for the entire face
2. Water container
3. Basic unit
4. Maximum filling level
5. Rotary switch
6. Indicator light
7. Aroma attachment
8. Facial sauna steam mask for nose and mouth
9. Steam valve (reverse)
10. Steam disc
11. Measuring cup

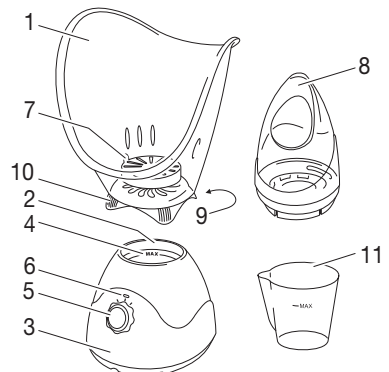6. Using the device

This device is intended only for use on the human face.

Using the facial sauna provides an ideal basis for cosmetic facial care. The skin of your face is prepared optimally for further treatments such as deep cleansing, moisturising, face masks etc.

Beauty salons very often begin facial treatments with a steam treatment. The invigorating steam increases the skin circulation, the pores are opened and impurities are removed more easily.

You can thus provide the perfect basis for every kind of facial care. The skincare products penetrate deep into the skin through the opened pores, thus optimising their action.

6.1 Operation

Place the device on a flat, firm, level and water-insensitive surface or on a table.

Use the measuring cup to fill the device's water container with water. Do not exceed the **maximum** filling level.

You should never fill the container directly under running water but only use the measuring cup that comes with it. The water volume at maximum filling level is enough for about 20 minutes use.

Do not use oils, herbs or other aromatic substances as a substitute for water but only added to adequate water filling.

As well as the facial mask, there is also a special attachment for the nose.

Fix the desired steam mask onto the basic unit by turning it.

Ensure that the mask engages firmly with the knobs.

6.2 Using the facial sauna

Uncoil the cable fully. Plug into the socket. Then turn the switch clockwise until the red indicator lights up. Turn the switch first to position "2" in order to speed up heating. After about 5 minutes, steam begins to rise. Turn the switch back to position "1" to reduce the steam output.

Adjustable steam outlet:

You can also regulate the strength of the steam outlet by turning the steam disc.

When the steam outlet is closed, the steam comes out of the steam valve which is located on the back of the device or mask. The water heating phase can also be shortened by closing the steam disc.

To avoid burns, never touch the steam outlet directly and ensure that the steam can come out of the valve freely.

Now sit down in front of the device. Bend forward until your face is above the mask.

If the steam effect is perceived as too hot, increase the distance to avoid burns.

Caution:

During the treatment, you should always turn the switch back to position "1".

Do not keep your face too close to the mask as otherwise burns can occur because of the intensity of the steam.

When you want to stop using the device or the water has evaporated, turn the switch back to position "0". The device's indicator goes out.

Rinse your face with cold water after use. Then apply a face cream or mask. The opened pores let the active ingredients penetrate deep into the skin.

Do not use the facial sauna for more than 20 minutes per day. Using it 2-3 times a week is recommended.

6.3 Use of medicinal herbs, oils and other aromatic substances

The aroma attachment enables you to use medicinal herbs, oils and other aromatic substances as part of an aromatherapy treatment. This is the best way for medicinal herbs, oils and other aromatic substances that mix with the steam to exert their effect. However, no solutions containing salt should be used.

Only use medicinal herbs, oils and other aromatic substances that are intended and suitable for use with a steam mask and are licensed for this by the producer of these substances.

Use:

Open the lid of the aroma attachment. Put a few drops of the oil on a cotton wool ball. Put in the cotton wool ball and then close the aroma attachment.

6.4 Using the nose/mouth steam mask

As the mucous membrane of the human nose is very sensitive, never use the device on position "2" when using with the nose/mouth steam mask. Do not keep your nose too close to the mask during use.

The treatment should not last more than 10 minutes per day and should be done 2–3 times a week at the most.

Use as described under 6.2.

7. Cleaning and care of the device

Unplug the device. Let it cool for about 10–15 minutes before taking the steam mask off the basic unit by turning it anticlockwise. Empty the water container by holding it upside down.

Clean grease and chalk deposits after each use with a commercially available decalcification agent suitable for aluminium.

Use a soft and slightly moistened cloth to clean the surface of the device and the steam mask and then wipe the surface dry.

- Do not use harsh agents (e.g. alcohol, acetone, benzine or abrasives) to clean the plastic parts.
- Unplug before cleaning.
- Do not use any sharp or pointed objects for cleaning.
- Do not immerse the device in water.

8. Technical specifications

Mains connection 230V ~ 50Hz 120W

9. Disposal

Please dispose of the device in accordance with the directive 2002/96/EG – WEEE (Waste Electrical and Electronic Equipment). If you have any queries, please refer to the local authorities responsible for waste disposal.

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