
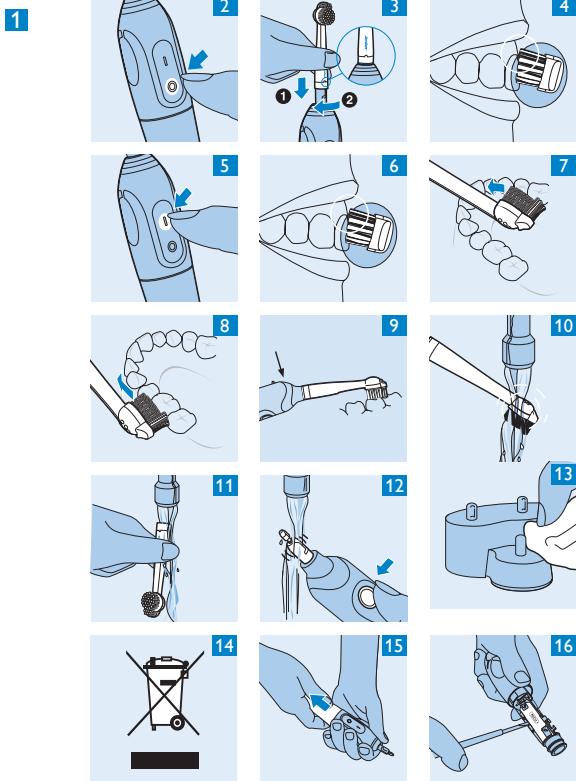
HX1630, HX1624, HX1622,
HX1620

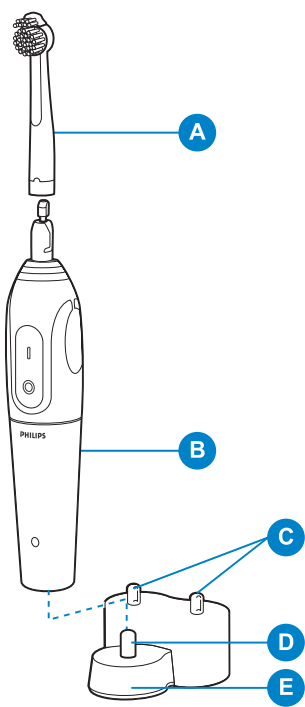PHILIPS

General description (Fig. 1)

- A Brush head
- B Handle
- C Brush head storage pins
- D Charger pin
- E Charger

Important

Read this user manual carefully before you use the appliance and save it for future reference.

- ▶ Check if the voltage indicated on the bottom of the charger corresponds to the local mains voltage before you connect the appliance.
- ▶ Keep the charger away from water. Do not place or store it over or near water contained in a bathtub, washbasin, sink etc. Do not immerse the charger in water or any other liquid.
- ▶ Never use the charger if it is damaged or if the cord or plug is damaged.
- ▶ The cord of the charger cannot be replaced. If the cord is damaged, the charger must be discarded.
- ▶ Make sure the socket is live when you recharge the toothbrush. The electricity supply to sockets in mirror cabinets over washbasins or shaver sockets in bathrooms may be cut off when the light is turned off.
- ▶ If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use the toothbrush.
- ▶ Do not use other brush heads than the ones recommended by the manufacturer.
- ▶ Replace the brush head every 3 months or sooner if signs of wear appear, such as crushed or bent bristles.
- ▶ This appliance is not intended for use by persons (including children) with reduced physical sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- ▶ Children should be supervised to ensure they do not play with the appliance.

Electromagnetic fields (EMF)

This Philips appliance complies with all standards regarding electromagnetic fields (EMF). If handled properly and according to the instructions in this user manual, the appliance is safe to use based on scientific evidence available today.

Preparing for use

Charging the appliance

- 1 Charge the appliance for at least 16 hours before you use it for the first time.
- ▶ Make sure the appliance is switched off when you start to charge it (Fig. 2).
- ▶ When the battery is being charged, the green indicator light is on (HX1630 only).
- ▶ You can leave the toothbrush in the charger and let it charge continuously.

Optimising the lifetime of the battery

After you have charged the toothbrush for the first time, do not recharge it between brushings. Continue to use it until the battery is empty before you recharge it again.

- Recharge the appliance when the motor almost stops running.
- Repeat this procedure three times in a row.
- To extend the lifetime of the battery, repeat this procedure every 6 months.

Using the appliance

- 1 To mount the brush head on the toothbrush, insert it into the bayonet catch and turn it clockwise a quarter turn (Fig. 3).

This prevents the brush head from becoming detached due to vibration.

- 2 Put toothpaste on the brush head.
- 3 Place the brush head in your mouth (Fig. 4).
- 4 Press the on/off button to switch on the toothbrush and start to brush your teeth (Fig. 5).

- ▶ For the best brushing method and the minimum brushing time, see 'Brushing tips'.

Brushing tips

Brushing method

- 1 Place the brush head against the tooth in such a way that it covers the surface of the tooth as well as the gum line (Fig. 6).
 - 2 Hold the brush head in this position for a few seconds and then move on to the next tooth (Fig. 7).
- The round brush head removes plaque from the teeth and beneath the gum line.
The individually moving Active Tip cleans between the teeth.
- 3 Guide the brush head carefully along the teeth and molars (Fig. 8).

Tip: Make sure you clean all surfaces of the teeth (inside, outside and chewing surface) as well as the gum line properly. A fixed brushing pattern prevents you from skipping certain areas.

Tip: For optimal dental hygiene, use dental floss and toothpicks in addition to the electric toothbrush.

Minimum brushing time

A brushing time of 2 minutes is generally long enough to clean your teeth properly. If you want, you can brush longer than 2 minutes for an even more thorough clean.

- ▶ A red indicator light will blink to indicate that 2 minutes have passed (HX1630 only).

Note: If plaque is removed properly, the gums will have a pink colour and will firmly surround the teeth. Healthy gums do not bleed when you eat or brush your teeth.

Note: If your gums bleed slightly the first few times you use this toothbrush, this bleeding will generally stop after a few days. If your gums continue to bleed, consult your dentist or dental hygienist.

Gum protection

- 1 The Gum Protection system guarantees that you will not exert too much pressure on the teeth and gums.
- 2 If you exert too much pressure on the toothbrush, you will feel and hear a click when the brush head flexes backwards. As soon as you reduce the pressure, the brush head will return to its normal position. This will protect your gums (Fig. 9)

Note: This will probably happen a few times when you first start using the brush. After some time you will have discovered the optimal brushing pressure and the brush head will flex backwards less often.

Cleaning and storage

- 1 Switch on the toothbrush and rinse the brush head under a hot tap (max. 60 °C). (Fig. 10)
- 2 Switch off the toothbrush. Detach the brush head and rinse the inside of the brush head (Fig. 11).
- 3 Switch on the toothbrush and rinse the drive shaft under the tap to remove toothpaste residues (Fig. 12).
- 4 Store the brush head on the brush head storage pin on the charger.
- 5 Clean the charger with a moist cloth (Fig. 13).

Unplug the charger before you start to clean it.

Replacement

- 1 Replace the brush head every 3 months for the best brushing results.
- ▶ Also replace the brush head when the bristles start to bend outwards.

You can obtain new brush heads from the shop where you purchased the toothbrush.
All round Philips replacement brush heads can be used.

Environment

- ▶ Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment (Fig. 14).
- ▶ The built-in rechargeable battery contains substances that may pollute the environment. Always remove the battery before you discard the appliance and hand it in at an official collection point. Dispose of the battery at an official collection point for batteries. If you have trouble removing the battery, you can also take the appliance to a Philips service centre. The staff of this centre will remove the battery for you and will dispose of it in an environmentally safe way.

Disposal of the battery

Make sure the battery is completely empty when you dispose of it.

- 1 Let the appliance run until it stops before you remove the battery.
- 2 Open the handle. (Fig. 15)
- 3 Use a screwdriver to pry the battery out of the battery holder (Fig. 16).

The toothbrush is no longer watertight once the housing of the toothbrush has been open. When you open the housing, your guarantee becomes invalid.

Guarantee & service

If you need information or if you have a problem, please visit the Philips website at www.philips.com or contact the Philips Customer Care Centre in your country (you find its phone number in the worldwide guarantee leaflet). If there is no Customer Care Centre in your country, go to your local Philips dealer or contact the Service Department of Philips Domestic Appliances and Personal Care BV.

Troubleshooting

If you are unable to solve your problem with the information below, please contact the nearest Philips service centre or the Philips Customer Care Centre in your country. You can also contact us online at www.philips.com.

Problem	Solution
The toothbrush does not work at all, although it has been in the charger overnight.	Make sure that the plug of the charger is properly inserted into the wall socket.
	Make sure the socket is live (see chapter 'Important').
	Make sure the toothbrush is switched off and recharge it for at least 16 hours.
The motor of the toothbrush works, but the brush head doesn't.	Check if the brush head has been properly attached to appliance (see also chapter 'Using the appliance').
The brushing power of the toothbrush has greatly decreased.	The toothbrush must be recharged. If you normally leave the toothbrush in the charger continuously, you may have to extend the lifetime of the battery as described in chapter 'Preparing for use'.
The toothbrush doesn't clean my teeth as well as before.	Replace the brush head.

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